

TACC2 Polyclonal Antibody

Purified Rabbit Polyclonal Antibody (Pab)
Catalog # AP58302

Specification

TACC2 Polyclonal Antibody - Product Information

Application IHC-P Primary Accession 095359

Reactivity Rat, Dog, Bovine

Host Rabbit
Clonality Polyclonal
Calculated MW 309427

TACC2 Polyclonal Antibody - Additional Information

Gene ID 10579

Other Names

Transforming acidic coiled-coil-containing protein 2, Anti-Zuai-1, AZU-1, TACC2

Format

0.01M TBS(pH7.4), 0.09% (W/V) sodium azide and 50% Glyce

Storage

Store at -20 °C for one year. Avoid repeated freeze/thaw cycles. When reconstituted in sterile pH 7.4 0.01M PBS or diluent of antibody the antibody is stable for at least two weeks at 2-4 °C.

TACC2 Polyclonal Antibody - Protein Information

Name TACC2

Function

Plays a role in the microtubule-dependent coupling of the nucleus and the centrosome. Involved in the processes that regulate centrosome-mediated interkinetic nuclear migration (INM) of neural progenitors (By similarity). May play a role in organizing centrosomal microtubules. May act as a tumor suppressor protein. May represent a tumor progression marker.

Cellular Location

Cytoplasm. Nucleus. Cytoplasm, cytoskeleton, microtubule organizing center, centrosome

Tissue Location

Strongly expressed in heart, skeletal muscle, brain, prostate, thyroid and trachea.

TACC2 Polyclonal Antibody - Protocols

Provided below are standard protocols that you may find useful for product applications.






• Western Blot

- Blocking Peptides
- Dot Blot
- Immunohistochemistry
- Immunofluorescence
- <u>Immunoprecipitation</u>
- Flow Cytomety
- Cell Culture

TACC2 Polyclonal Antibody - Images

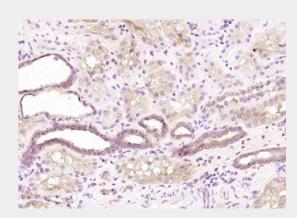

Paraformaldehyde-fixed, paraffin embedded (Human kidney); Antigen retrieval by boiling in sodium citrate buffer (pH6.0) for 15min; Block endogenous peroxidase by 3% hydrogen peroxide for 20 minutes; Blocking buffer (normal goat serum) at 37°C for 30min; Antibody incubation with (TACC2) Polyclonal Antibody, Unconjugated (bs-5736R) at 1:200 overnight at 4°C, followed by operating according to SP Kit(Rabbit) (sp-0023) instructions and DAB staining.