

Anti-Rabbit IgG (H&L) (Peroxidase Conjugated) Secondary Antibody

Chicken Polyclonal, Peroxidase (Horseradish)
Catalog # ASR2490

Specification

Physical State

Anti-Rabbit IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Product Information

Description Anti-RABBIT IgG (H&L) (CHICKEN) Antibody

Peroxidase Conjugated

Host Chicken

Conjugate Peroxidase (Horseradish)

Target Species Rabbit
Reactivity Rabbit
Clonality Polyclonal
Application ,1,10,15,

Application Note ELISA 1:140,000; Western Blot

1:2,000-1:10,000;Immunochemistry

1:1,000-1:5,000 Lyophilized

Host Isotype IgG

Target Isotype IgG (H&L)

Buffer 0.02 M Potassium Phosphate, 0.15 M

Sodium Chloride, pH 7.2

Immunogen Anti-Rabbit IgG whole molecule was

produced by repeated immunization with Rabbit IgG whole molecule in chicken.

Reconstitution Volume 1.0 mL

Reconstitution Buffer Restore with deionized water (or

equivalent)

Stabilizer 10 mg/mL Bovine Serum Albumin (BSA) -

Immunoglobulin and Protease free

Preservative 0.01% (w/v) Gentamicin Sulfate. Do NOT

add Sodium Azide!

Anti-Rabbit IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Additional Information

Shipping Condition

Ambient

Purity

This product was prepared from monospecific antiserum by immunoaffinity chromatography using Rabbit IgG coupled to agarose. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Peroxidase, anti-Chicken Serum, Rabbit IgG and Rabbit Serum.

Storage Condition

Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.



Precautions Note

This product is for research use only and is not intended for therapeutic or diagnostic applications.


Anti-Rabbit IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Protein Information

Anti-Rabbit IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Protocols

Provided below are standard protocols that you may find useful for product applications.

- Western Blot
- Blocking Peptides
- Dot Blot
- Immunohistochemistry
- Immunofluorescence
- <u>Immunoprecipitation</u>
- Flow Cytomety
- Cell Culture

Anti-Rabbit IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Images

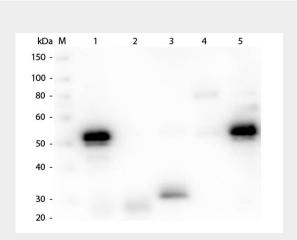

Western Blot of Anti-Rabbit IgG (H&L) (CHICKEN) Antibody . Lane M: 3 μ l Molecular Ladder. Lane 1: Rabbit IgG whole molecule . Lane 2: Rabbit IgG F(ab) Fragment . Lane 3: Rabbit IgG F(c) Fragment . Lane 4: Rabbit IgM Whole Molecule . Lane 5: Normal Rabbit Serum . All samples were reduced. Load: 50 ng per lane. Block: MB-070 for 30 min at RT. Primary Antibody: Anti-Rabbit IgG (H&L) (CHICKEN) Antibody 1:1,000 for 60 min at RT. Secondary antibody: Anti-Chicken IgG (GOAT) Peroxidase Conjugated Antibody 1:40,000 in MB-070 for 30 min at RT. Predicted/Obsevered Size: 25 and 50 kDa for Rabbit IgG and Serum, 25 kDa for F(c) and F(ab), 70 and 23 kDa for IgM. Rabbit F(c) migrates slightly higher.

Anti-Rabbit IgG (H&L) (Peroxidase Conjugated) Secondary Antibody - Background

Anti-Rabbit IgG whole molecule antibody generated in chicken detects specifically Rabbit IgG whole molecule. This secondary antibody anti-Rabbit is ideal for investigators who routinely perform ELISA, Sandwich ELISA, titration assays, western-blot, immunoprecipitation and more generally immunoassays.